

**A complete solution to minimallyinvasive
disposable hysteroscopy**



• Office Solution



HTx40 is compatible with all 5Fr rigid/semirigid instruments without using speculums.

• Operating Room Solution



HTx60/60s/60T is compatible with 7Fr and 9Fr instruments including 9Fr tissue removing devices.



Automatic white
balance/gain/focus




Integrated processor and display
with a built-in battery enabling
portability and versatility




on-board patient
management system with
footage transfer capabilities



13.3", 1080p medical display
monitor with patented digital
image processing

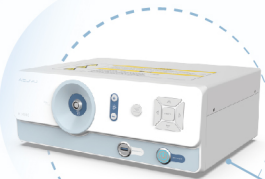

HTx1000

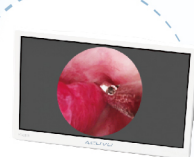

HTx2000



AcuVu's aims to advance gynecology by offering a disposable hysteroscope platform that enhances patient comfort, safety, and outcomes. Our mission is to provide accessible and innovative solutions that empower healthcare providers to deliver top-notch care in any office setting.

MORE INFO AT:

www.acuvuinc.com

**SEE CLEARLY
TREAT SIMPLY**



HummingScope®

**DISPOSABLE DIAGNOSTIC & OPERATIVE
HYSTEROSCOPY**





SEE CLEARLY TREAT SIMPLY



No capital investment



See and treat in office



No risk of cross-contamination



4.5mm OD, No anesthesia
and speculum needed



Compatible with all 5Fr, 7Fr,
9Fr flexible/rigid/semi-rigid instruments




Intuitive orientation in uterus, excellent fluid circulation,
and effortless photo/video recordings




MODEL	HTx 40	HTx 60	HTx 60s	HTx 60T
Outer Diameter	4.5mm	6.2mm	6.2mm	6.2mm
Working Channel	2.1mm	3.3mm	3.3mm	3.3mm
Instruments Compatibility	5Fr flexible/ semi-rigid	7Fr, 9Fr rigid	MyoSure Lite/ Reach/Manual	TruClear mini
Working Length	240mm	240mm	195mm	240mm
Field of View	120 °	120°	120°	120°
Direction of View	8°	12°	0°	0°



Tightly sealed
without any leakage
with tools

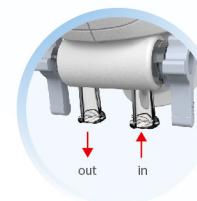

Separate in and out water channels
for better fluid management

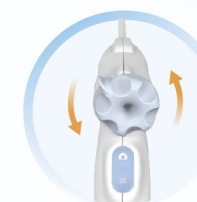

360 degrees knob rotation for
complete uterus visibility



Exceptional image quality
that fulfills both diagnostic
and therapeutic purposes