

# Litavo


**09221** The ideal aid while reading or eating a meal, both in bed and chair. The tabletop can be adjusted continuously and is equipped with a tiltable part that can be used to read on.

## Base


- Firmly welded U-shaped frame.
- The high/low column consists of 2 aluminium profiles and can be telescopically adjusted in height by means of a gas spring.
- The column can be installed at the left or right side, depending on the customer's wishes.
- The gas spring is activated by means of a lever near the fixed tabletop and guarantees a continuous height adjustment of circa 69 cm to 99 cm.
- The U-shaped base can be optionally equipped with two stainless steel protection caps.
- The metal support is protected with an epoxy-coating, *RAL 9006 blank alu.*
- The aluminium parts are protected by an anodised layer.
- Equipped with 4 design castors Ø 65 mm, two of which have brakes.

## Tabletop

- The tabletop consists of a fixed and a tiltable part. Both parts are equipped with an upright synthetic edge.
- The tiltable part is supported by an aluminium profile.
- The tabletop can be tilted in 2 positions: 22,5° and 45°.
- The tabletop is tilted by means of an easily accessible lever on the outside.
- The tabletop is made of solid laminate, thickness 8 mm.

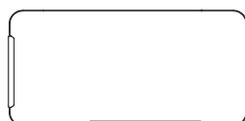


One-piece tabletop symmetrical version

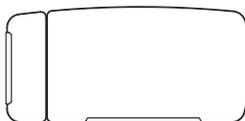

Two-piece tabletop left version



Two-piece tabletop right version


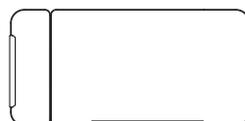


Two-piece tabletop symmetrical version



## Dimensions

**09221** | 89 x 43 x H 69-98 cm  
Gap distance U-shape base: 68 cm  
Height base: 10 cm

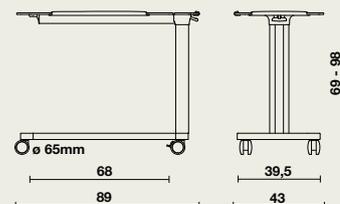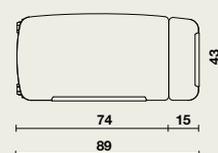

## Materials and finishing

- Metal: epoxy coating RAL 9006 blanc alu.
- Solid laminate.
- Synthetic.

For further information about used materials, constructions and maintenance: consult our materialfiles.

## Options

- Tabletop consists of a one piece fixed part (non tiltable, symmetrical).
- Stainless steel protection caps.
- 2 suspension hooks.
- Larger base of 12 cm, gap distance: 80 cm.
- Symmetrical table top.
- Gallery on fixed part (only possible with symmetrical tabletop).
- Reading support on the second long side of the tabletop (only possible with symmetrical tabletop).
- Reading support on the second short side of the tabletop (only possible with symmetrical tabletop).

## Colours

We would like to refer to our extensive colour gamut.  
Sample card + samples available on request.

