

Reading table type Livia

04202 Livia is the ideal patient aid to read or eat a meal while lying in bed or sitting in a chair. The tabletop can be continuously adjusted and is equipped with a folding part that can be used as a reading stand.

Base

- Firmly welded U-shaped frame.
- The high/low column is supported by means of a gas spring.
- The gas spring is activated by means of a pedal and is continuously adjustable between 82 cm and 115 cm.
- Depending on the customer's wish, the column can be installed at the right or left side of the base.
- The base can be optionally equipped with two stainless steel cover caps.
- All parts of the tabletop mechanism are made of high-quality synthetic materials or aluminium.
- All metal parts are protected by an epoxy coating.
- The aluminium guide profile is equipped with an anodised layer.
- 4 design castors Ø 65 mm, 2 of which have brakes.

Tabletop


- The tray is made of fire-retardant ABS.
- Equipped with an upright edge.
- The tray can be vertically stowed.
- When folded, the tabletop can be easily inclined.
- The inclination is possible in both directions, and in 2 positions: 22,5° and 45°.



Dimensions

04202 | Reading table type Livia

- 70 x 40 x 82-115 cm
- Tabletop: 63 x 40 cm

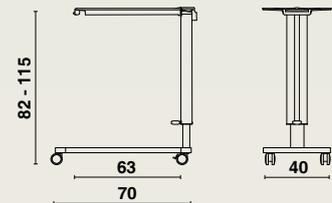

Materials and finishing

- Metal: Steel, epoxy coating RAL 9006, anodised aluminium.
- Synthetic.

For further information about used materials, constructions and maintenance: consult our materialfiles.

Options

- Stainless steel cover caps.

