

# Atlas Care

**06261** The Atlas Care bedside table with refrigerator, is very compact and matches perfectly with the Atlas wardrobe gamut.

## Version with refrigerator (refrigerator, open compartment, drawer)

- 2 side panels with a front lath, back, bottom, and a tabletop.
- Made of melaminated chipboard (thickness 18 mm).
- The panels are connected to each other by means of wooden dowels and eccentric cams.
- Storage cabinet with refrigerator.
- Noiseless absorption refrigerator "Dometic Minicool - Silencio DS 400 BI".
- Capacity: 37 litres.
- Door with HPL finishing.
- Interior lighting with LED.
- Spiral cable.
- The removable drawer is made of metal sides, between which a wooden bottom and back have been attached.
- The drawer is installed in the cupboard frame by means of high-quality steel rails with an end stop.
- Metal design handle.
- The metal parts are protected by an epoxy-coating.
- The tabletop has rounded corners and edges and is made of melaminated chipboard of 18 mm.
- 4 design twin castors Ø 65 mm.


## Optional dinner tray Livia

- Folding.
- Easy to incline in both directions.
- The dinner tray has an upright edge.
- The folding part is fixed on a strong aluminium guiding profile that has been finished with an anodisation layer.
- The dinner tray is height-adjustable (between 80 and 110 cm) by means of a gas spring.
- Dimensions of the dinner tray: 63 x 40,5 cm.
- The gas spring is activated by means of a foot pedal.
- The dinner tray is made of fire-retardant ABS.
- The bedside table stands on 5 design castors Ø 65 mm.
- The dinner tray can be optionally provided with a Click 'n' turn mechanism which allows to move the dinner tray from one side to the other side of the bedside table.



## Dimension

**06261** | 49,5 x 49,5 x height 88 cm  
width dinner tray Livia | 65 x 49,5 x height 88 cm

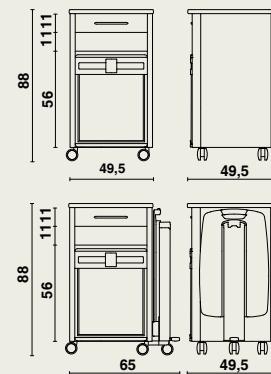

## Materials and finishing

- Wood: chipboard (melaminated).
- Metal: epoxy-coating (clear gray).
- Withstands the common cleaning products.

For more information about the **materials, construction and maintenance**: please consult our material lists.

## Options

- Finishing in veneer.
- Livia dinner tray.
- Gallery.
- Cylinder lock
- Cases for accessories on the tabletop.
- Telephone holder.
- Extensible tray.
- 4 design castors Ø 65 mm.
- 4 design castors Ø 65 mm with brake.
- Click 'n' turn mechanism.
- USB charger and 2 sockets built into the back.

## Colours

We would like to refer to our extensive colour gamut.  
Sample card + samples available on request.

## Handles

- Open handle
- Closed handle

