





WEL Communications

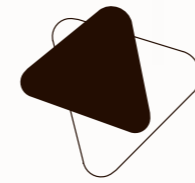

# Siena & Aris Care

by Haelvoet

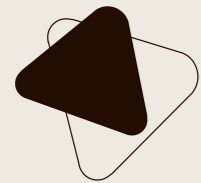

HAEVOET NV **BELGIUM**  
L. BEKAERTSTRAAT 8, 8770 INGELMUNSTER  
T +32 (0) 51 48 66 95 F +32 (0) 51 48 73 19  
INFO@HAEVOET.BE WWW.HAEVOET.BE

HOSPITAL - NURSING HOME FURNITURE

# Siena

The Siena nursing bed combines user-friendliness with a highly functional and efficient mattress platform.




001	<b>Sideguards (*) • Safety has no price</b>
	The Atmosphère sideguards protect the patient over 2/3 of the bed length, whereas the Pollux or Trix sideguards offer maximum patient protection over the entire bed length.
002	<b>Mattress platform • For a comfortable lying position</b>
	You can choose between 2 different mattress platforms. The standard version is equipped with well-ventilating metal mattress platform profiles. The optional solid laminate mattress platform sections (*) can be easily removed and cleaned.
003	<b>Bed panels • A harmonious duo</b>
	You can choose between 4 different bed panel types. All types give the bed a tasteful wooden touch.
004	<b>Mattress holders • Purposefully present</b>
	The mattress holders prevent the mattress from moving towards the bed sides. They can also be used as a handle to adjust the backrest and knee rest manually.
005	<b>Cases • Optimal care</b>
	Two cases are provided at the head of the bed, and can be used to install an IV rod, a lifting pole or other accessories.
006	<b>Operating lever • Central and within reach</b>
	The operating lever activates a gas spring, so that the desired mattress platform section can be adjusted. Optionally, the backrest and legrest adjustments can be activated electrically (*). It is also possible to switch off the electric functions with one single button on the lock device (*).
007	<b>Accessory holders (*) • First-class functionality</b>
	Three multifunctional holders on both sides make it possible to install accessories very fast and easily.
008	<b>Wooden surrounds of the mattress platform frame (*) • An extra trump</b>
	The wooden surrounds on the mattress platform frame give the bed a homely look.
009	<b>Bed extension • Increased flexibility</b>
	The standard integrated bed extension makes it possible to extend the bed with 20 cm. The optional footrest extension (*) guarantees a complete functionality of the mattress platform.
010	<b>Castors (*) • Superb manoeuvrability</b>
	Our extensive gamut of high-quality design castors guarantees a highly mobile and manoeuvrable bed.

## Comfort for everyone

Nowadays, both the resident and nursing staff demand comfortable and ergonomic beds. Only then, the resident can be as autonomous as possible, whereas the nursing staff can work efficiently.



### • Activations

The operating lever is user-friendly and easily accessible on both sides of the bed.

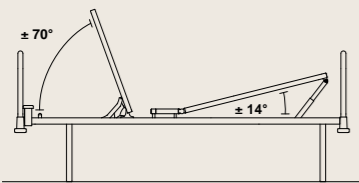

#### Seating or relax position

The three- or four-section (\*) mattress platform offers ample adjustment possibilities and can be easily activated on both sides of the bed. Thanks to this, the desired seating position – or relax position (\*) in case of a four-section mattress platform – can be easily obtained. Optionally, the backrest and legrest adjustments can be activated electrically by means of a hand switch (\*). You can also opt for a lock device (\*) which makes it possible to switch off several electric bed functions with one single button.

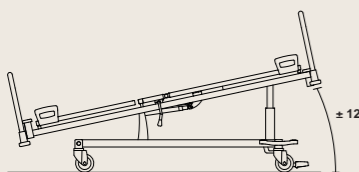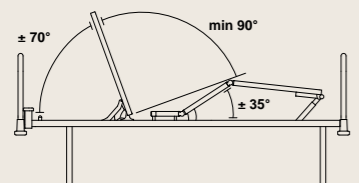


#### Trendelenburg position

The Trendelenburg position (\*) can be activated by means of the red pedal at the foot of the bed. When the pedal is stepped on, the bed inclines fluently until the desired position is obtained.

### • Retracting backrest

A retracting backrest increases the patient's lying comfort considerably. While inclining the backrest upwards, it is automatically moved backwards. By doing so, pressure is relieved on the pelvis, whereas the bedsore risk is reduced substantially.

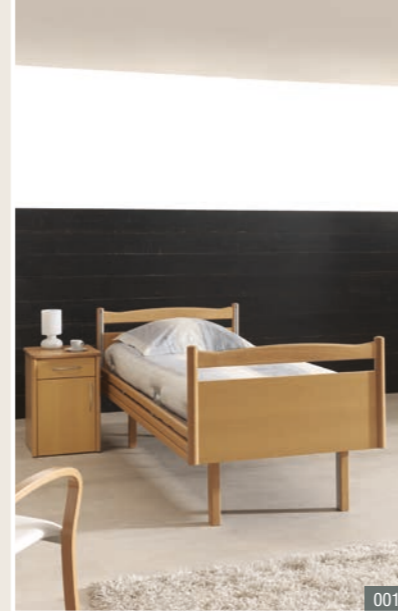

### • High/low carriage

The high/low carriage Vico (\*) can be installed very easily to obtain a complete high/low bed. Combined with the removable synthetic cover cap (\*), even the surfaces that are difficult to reach can be cleaned easily.

## User-friendly

Particular attention has been paid to efficient care during the development of the Siena bed.

### • Castors - Excellent mobility

The daily use of the Siena bed is simplified, thanks to its excellent mobility and manoeuvrability. You can choose from our extensive gamut of castors (Ø 100 mm or Ø 125 mm) with or without central brakes (\*).



standard on fixed legs or height-adjustable (\*) legs



castors Ø 100 mm with individual brakes, optionally height-adjustable (\*)



castors Ø 125 mm with central brake and directional castor (\*)  
The directional castor makes it easy to move the bed in a straight line.




### • Bed extension

A nursing bed has to meet everyone's needs, irrespective of the patient's length. The integrated bed extension makes it possible to extend the bed with 20 cm.

Two different methods are possible to extend the mattress platform. You can choose between telescopic footrest extension combined with the extendible Trix sideguards or a removable footrest extension of 20 cm.

## Safety galore

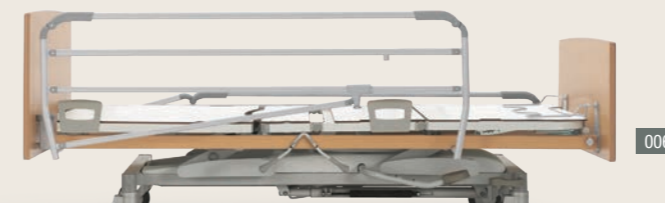

001 The Siena Care bed combined with the Atlas bedside table.

002 The operating lever.

003 The Trendelenburg pedal (\*).

004 The retracting backrest of the mattress platform.

005 The Trix sideguards (\*) protect the resident over the entire bed length, thanks to the unique telescopic extension.

006 The Siena Care bed with solid laminate mattress platform sections (\*), Atmosphère sideguards (\*) and wooden surrounds (\*).


007 The Siena Care bed with steel mattress platform profiles, Pollux sideguards (\*), and wooden surrounds (\*).

008 The Siena Care bed with steel mattress platform profiles, Trix sideguards (\*), and wooden surrounds (\*).

\* option for Siena Care

# Siena Care Summary

## • Dimensions and performances (standard bed)

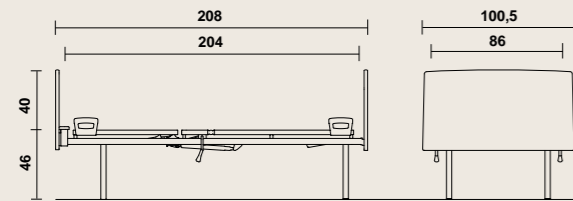

- length: 208 cm, width: 100,5 cm;
- height of the mattress platform: ± 46 cm, optionally adjustable up to ± 55 cm (in 3 positions);
- safe working load: 180 kg;
- weight (without accessories): 65 kg;
- mattress platform dimensions: 204 x 86 cm;
- mattress dimensions: 195 x 85 cm, thickness: 12 to 16 cm (\*).

\* The mattress thickness depends on the sideguards type.



Haelvoet makes sure that its excellent price-quality ratio is not at the expense of the environment. Our company considers durability – in all possible meanings – of paramount importance. All our products have a long lifespan and are delivered with a 10-year-guarantee on construction errors.



Choosing Haelvoet is choosing an honest product, thanks to our environment-friendly production processes, a maximum use of recyclable materials and our PEFC certificate.



The Siena nursing bed is made in accordance with ISO 9001:2015 and ISO 13485:2016 and complies with the European Medical Device Regulation EU 2017/745 and all other applicable European harmonised norms:

- EN 60601-1 Medical electrical equipment - Part 1: General requirements for basic safety and essential performance.
- EN 60601-1-2 Medical electrical equipment - Part 1-2: General requirements for safety and essential performance - Secondary norm: Electromagnetic compatibility - Requirements and tests.
- EN 60601-2-52 Medical electrical equipment - Part 2-52: Particular requirements for the basic safety and essential performance of medical beds, except for the electrical clauses.
- EN ISO 14971 Application of risk management to medical devices.


	Standard Siena		Options Siena				
MATRESS PLATFORM							
	• 3-section steel mattress platform	• Telescopic bed extension up to 20 cm	• 4-section steel mattress platform	• 3-section mattress platform with removable solid laminate sections	• 4-section mattress platform with removable solid laminate sections	• Footrest extension up to 20 cm	
CASTORS							
	• Legs, fixed height		• Adjustable legs in 3 positions	• Non-adjustable legs, castors Ø 100 mm, individual brakes	• Adjustable legs, castors Ø 100 mm, individual brakes	• Non-adjustable legs, castors Ø 125 mm, central brakes	• Non-adjustable legs, castors Ø 125 mm, directional castor, central brakes
PANELTYPE							
	• Midea panel		• Puro panel	• Felina panel	• Valentino panel	• Bumpers (not on Midea panel)	• Wall spacer
SIDEGUARDS							
			• Atmosphère sideguards	• Pollux sideguards	• Trix sideguards	• Trix sideguards with 4th siderail	
CONTROLS							
	• Operating lever		• Electric hand switch	• Lock device (ACL)	• Trendelenburg pedal (fixed height)		
EXTRAS							
	• 2 cases at the head end		• High/low carriage Vico	• Wooden surrounds of the mattress platform	• Wooden surrounds of the legs	• 6 synthetic accessory hooks	



# Aris

The Aris Care bed is an efficient and comfortable nursing bed. The optional inclining backrest can be easily operated on both bed sides.


001	<b>Mattress platform</b> • For a comfortable lying position
	The mattress platform is standard equipped with well-ventilating metal mattress platform profiles.
002	<b>Bed panels</b> • Stylish and durable
	The bed panels have a stylish look. A combination of a high head panel and a low foot panel is possible, as well as the combination of two high or two low panels.
003	<b>Cases</b> • Optimal care (only for a 2-section mattress platform)
	Two cases at the head can be used to install an IV rod, a lifting pole or other accessories.
004	<b>Operating lever (*)</b> • Central and within reach (only for a 2-section mattress platform)
	The operating lever on both sides activates a gas spring, which makes it easy to adjust the backrest in the desired position.
005	<b>Wooden surrounds of the mattress platform frame (*)</b> • An extra trump
	The wooden surrounds on the mattress platform frame give the bed a homely look.

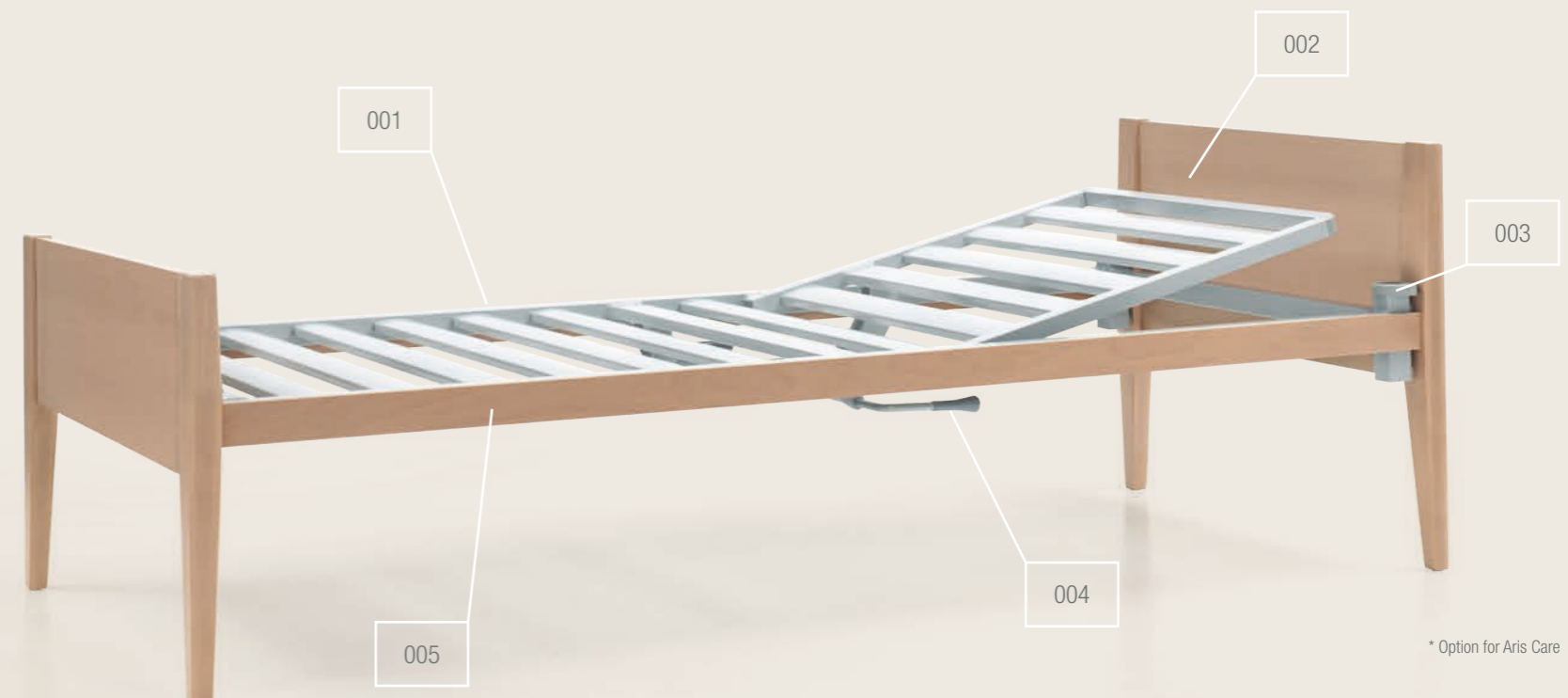

\* Option for Aris Care

# Aris Care Summary

## • Dimensions

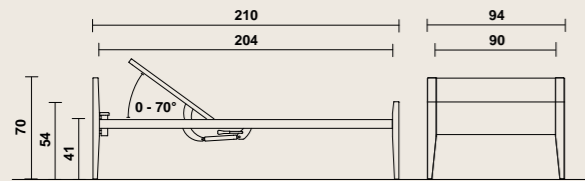

- length: 210 cm, width: 94 cm;
- height of the mattress platform: ± 41 cm;
- mattress platform dimensions: 204 x 90 cm;
- mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm.



Haelvoet makes sure that its excellent price-quality ratio is not at the expense of the environment. Our company considers durability – in all possible meanings – of paramount importance. All our products have a long lifespan and are delivered with a 10-year-guarantee on construction errors.



Choosing Haelvoet is choosing an honest product, thanks to our environment-friendly production processes, a maximum use of recyclable materials and our PEFC certificate.



The Aris nursing bed is made in accordance with ISO 9001:2015 and ISO 13485:2016 and complies with the European Medical Device Regulation EU 2017/745 and all other applicable European harmonised norms:

EN 60601-1 Medical electrical equipment - Part 1: General requirements for safety, except for the electrical clauses.

EN 60601-2-52 Medical electrical equipment - Part 2-52: Particular requirements for the basic safety and essential performance of medical beds, except for the electrical clauses.

EN ISO 14971 Application of risk management to medical devices.

	Standard Aris	Options Aris
MATTRESS PLATFORM	<ul style="list-style-type: none"> <li>• 1-section steel mattress platform</li> </ul>	<ul style="list-style-type: none"> <li>• 2-section steel mattress platform</li> </ul>
CASTORS	<ul style="list-style-type: none"> <li>• Legs with synthetic glides</li> </ul>	
PANELTYPE	<ul style="list-style-type: none"> <li>• Bed panel Aris Care</li> </ul>	<ul style="list-style-type: none"> <li>• Combination of high head panel/low foot panel</li> </ul>
CONTROLS		<ul style="list-style-type: none"> <li>• Operating lever</li> </ul>
EXTRAS	<ul style="list-style-type: none"> <li>• 2 cases at the head of the bed (only with a 2-section mattress platform)</li> </ul>	<ul style="list-style-type: none"> <li>• Wooden surrounds of the mattress platform</li> </ul>

# Accessories

The accessories we offer can be harmoniously connected with the existing cases and accessory holders. Haelvoet strongly pays attention to this, so that the accessories offer you the same top quality as our nursing beds. After all, only the best quality is good enough for our customers.

001	Lifting pole	<p>The lifting pole is made of a bent thick-walled steel tube. It is equipped with a handle and a ribbon that can be manually adjusted. An automatic roll-up system is optionally available as well.</p>	006	IV rod	<p>The IV rod consists of a fixed and a height-adjustable part that is equipped with 4 hooks. The IV rod is also available in a non-adjustable version with 2 hooks.</p>
002	IV rod for lifting pole	<p>The IV rod for the lifting pole consists of 4 hooks and can be easily installed on the lifting pole.</p>	007	Raising aid (*)	<p>The ergonomically-designed raising aid offers optimal patient support while getting in and out of the bed.</p>
003	Dinner tray for sideguards (*)	<p>The wooden dinner tray can be installed on all types of sideguards. It can be easily hung onto the foot panel.</p>	008	Urine bottle holder (*)	<p>The urine bottle holder is installed in two optional synthetic accessory hooks. You can choose between two different models.</p>
004	Accessory bracket (*)	<p>The accessory bracket can be used in case of specific customer needs. The bracket has 4 hooks and has a length of 310 mm. It is put in the optional accessory hooks.</p>	009	Cover for sideguards (*)	<p>The removable cover is made of imitation leather and can be easily put around the sideguards. This cover is especially recommended to protect restless residents.</p>
005	Bowl holder (*)	<p>Two bowls can be placed into the bowl holder. The holder is installed in the optional accessory hooks.</p>	010	Removable footrest extension (*)	<p>To use the mattress platform optimally when the bed is extended, you can opt for a removable footrest extension of 20 cm that has to be put at the foot end.</p>

\* only possible for Siena Care