


INFRARED Thera duo

LiD
Eclairage Lighting 

A two headed heating lamp with two 250 W bulbs, on two 50 cm long goosenecks that are very durable and silent. A session can be performed in half the time on two different zones (neck/shoulders - neck/back - the two thighs, etc.) with better heat distribution than with one lamp, or one treatment zone can be treated more intensely (500 W) with careful placing of the two heads.

Can also be used to act on all sides of a muscle group (e.g.: Biceps femoris, rectus femoris, vastus lateralis...).
Needing only a single plug, the Thera Duo is functionally mobile with a small footprint.

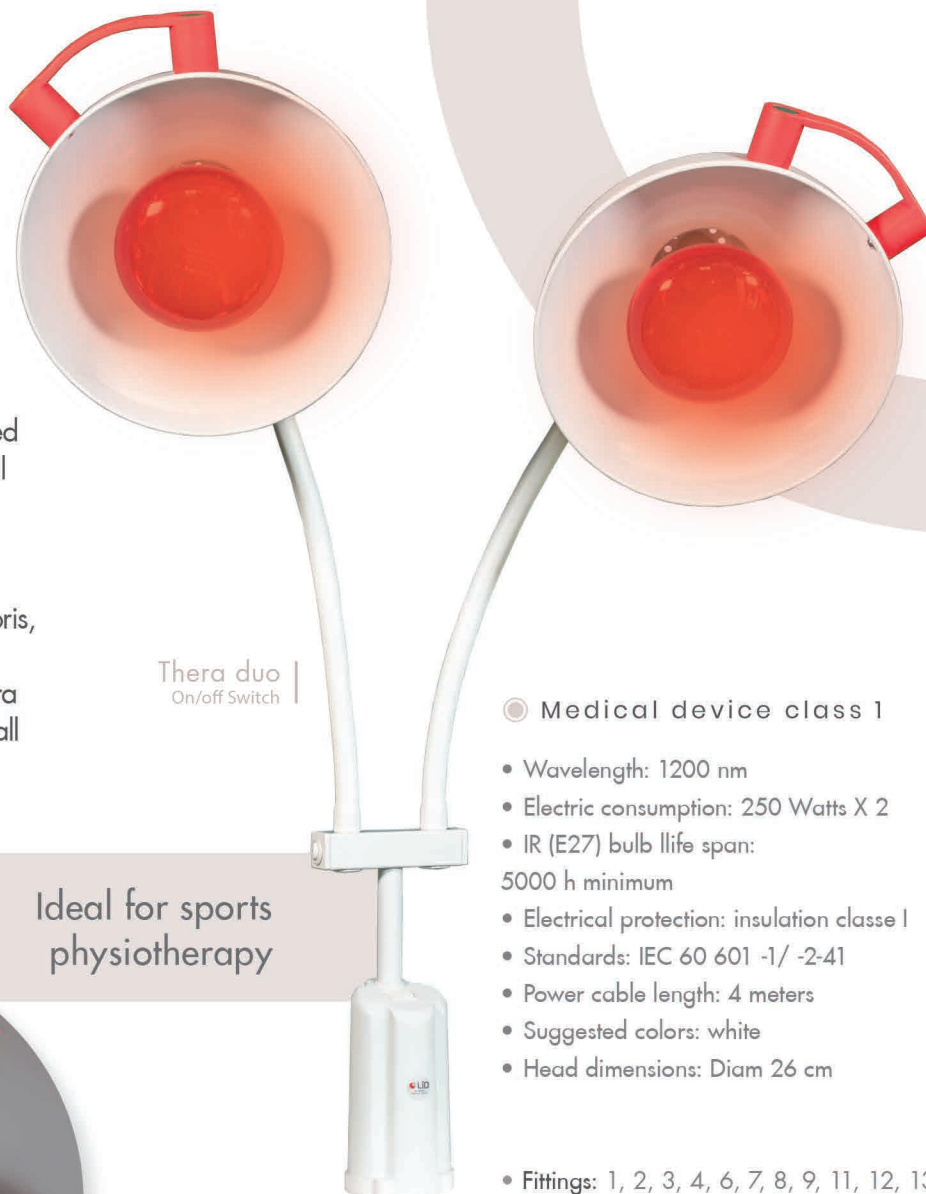

Thera duo |
On/off Switch

Ideal for sports
physiotherapy

● Medical device class 1

- Wavelength: 1200 nm
- Electric consumption: 250 Watts X 2
- IR (E27) bulb life span: 5000 h minimum
- Electrical protection: insulation classe I
- Standards: IEC 60 601 -1/ -2-41
- Power cable length: 4 meters
- Suggested colors: white
- Head dimensions: Diam 26 cm

- Fittings: 1, 2, 3, 4, 6, 7, 8, 9, 11, 12, 13



CODES

	On/off Switch	Contactless switch	Timer switch
Thera duo	IRTDUO500	IRTDUO500-ISC	IRTDUO500MI

IR bulb **87115005**



HIGH
HEAT OUTPUT