




Digital Video

Laryngoscope

touching for life

SmartScope VL Pro

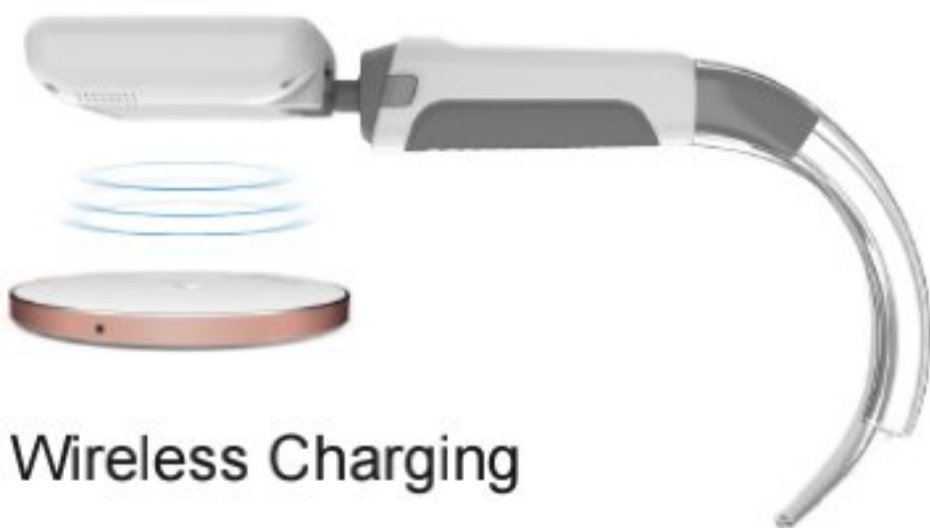





Wireless Charging

-  Safety
-  Fast
-  HD
-  Easy
-  Cost Effective



Five solution for the Intubation

 Safety	 Fast	HD High- Definition	 Easy	 Cost Effective
<ul style="list-style-type: none"> • Ergonomic and durable design for long time usage • High quality material with anti shatter protection for whole unit • Disposable Blades avoid cross infection • Anti- Fog Lens and Anti-Microbial Handle 	<ul style="list-style-type: none"> • Fast intubation after short time training • Fast install or remove disposable blade • Lightweight and portable 	<ul style="list-style-type: none"> • 3.5 inch anti- reflective full touch display with 640*480 resolution • 2.0 Mega Pixel Camera provide 1280*720P photo and video • Rotation angle for up-down $\geq 120^\circ$ and left-right $\geq 270^\circ$ 	<ul style="list-style-type: none"> • One shot button for Recording • One cable connect to PC for data management • One cable connect to large LCD for better view (option) • Full touch scree with user friendly operation • Wireless charging 	<ul style="list-style-type: none"> • Favorable cost for all kinds of disposable blades • Rechargeable battery for at least 4 hours working time

Suitable for

• Difficult incubation of glottic opening, easily guide the medical staff to accurately perform endotracheal intubation for anesthesia or first aid.



Anesthesia center



Respiratory department



ICU ward



Emergency center

Technical Specifications

Display	Monitor: LCD
	Resolution Ratio: 640X480RGB
	Aspect Ratio: 4:3
	Video refresh Rata: 30FPS
Data output	Easy to establish and store files in JPEG or MP4. 16GB standard can be upgrade to 256GB. Connect to PC by USB-Type C cable.
Image output	Connect to external LCD by HDMI cable, convenient for training and presentation
Camera and image	Resolution Ratio 2.0M Pixels
	Field Angle: $\geq 60^\circ$
	Illuminance $\geq 1000\text{LUX}$
	Working distance: 20mm-100mm
	Color temperature $\geq 5000\text{K}$
Power and Battery	Type Rechargeable Lithium Battery
	Battery life Cycle: More Than 500 times
	AC100-240V, 50/60Hz, 50VA
	Charging times: less than 4hrs
	5W Wireless Charging (Optional)
Working Condition	Temperature: $-10\sim+45^\circ\text{C}$
Size and weight	Product Size: 245*100*130mm; Product weight: 0.28kg;
	Carton size: 345*120*282mm; Carton weight: 2.2kg; Packing quantity: 1pcs



S4 Adult large(6.0mm or more)



S3 Adult (6.0~9.0mm)



S2 Child (3.0~6.0mm)

Specification subject to be changed without prior notice.